

---

# Perio Voice

---

## Installation and User Guide

---

Planet DDS


---

## Table of Contents


Installation .....	3
Requirements.....	7
Internet Explorer.....	7
Firefox .....	8
Chrome .....	10
User Guide .....	11
How to Configure.....	11
How to Run .....	12
How to Use .....	14
New Exam .....	14
Pocket Depth.....	14
FGM.....	14
MGJ .....	14
Bleeding .....	14
Suppuration.....	14
Furcation .....	14
Mobility .....	14
Next.....	14
Previous.....	14
Jump to #Tooth No .....	14
Go lower .....	14
Go upper .....	14
Voice Feedback .....	15
New updated version install Instructions .....	16

# Installation


1. Download **Denticon Voice Command** app from **Denticon** (Help > [Help & Support](#)).

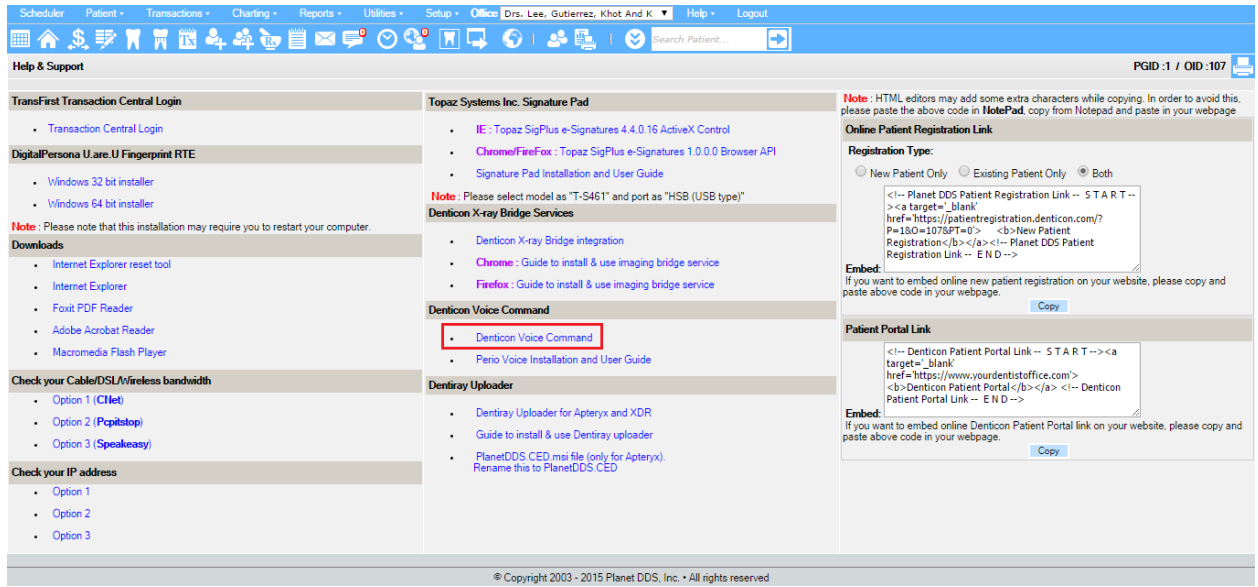


2. Right click on the downloaded file and click **Install**.

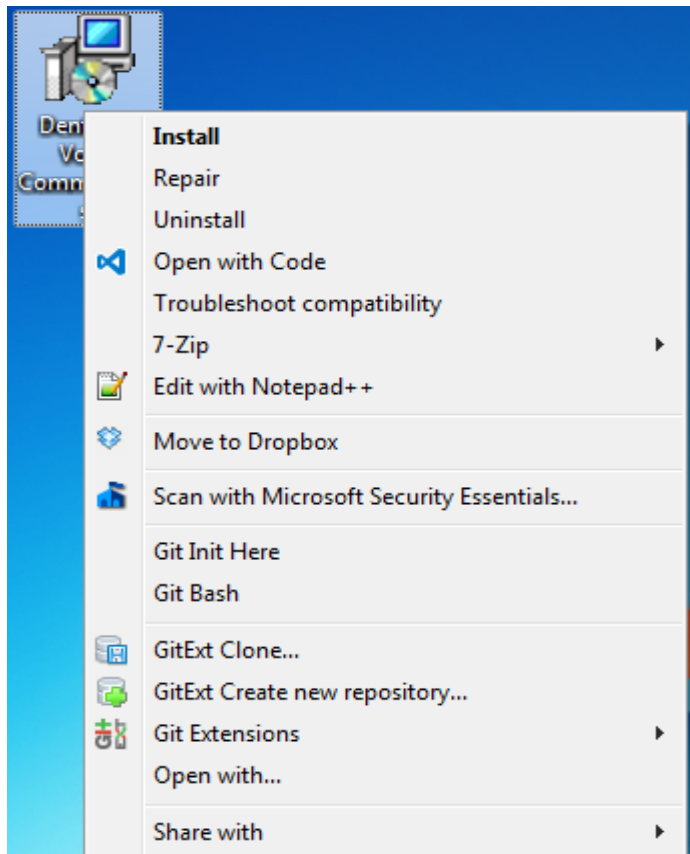


3. The installation will start. Click **Run**.

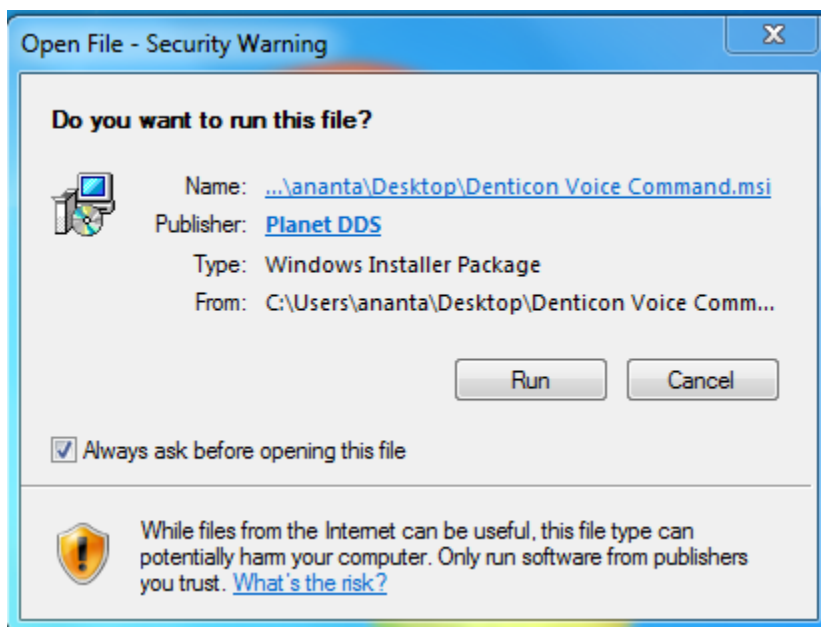


4. Click **Next**.

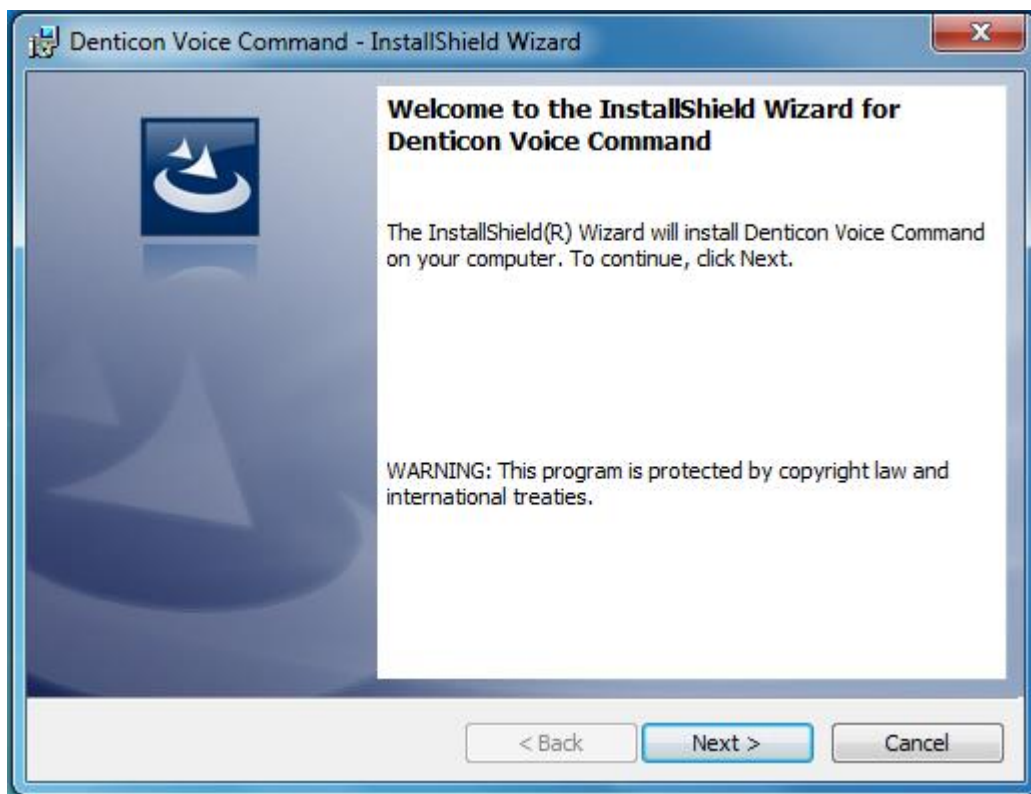


5. Select **"I accept the terms in license agreement"** and click **Next**.

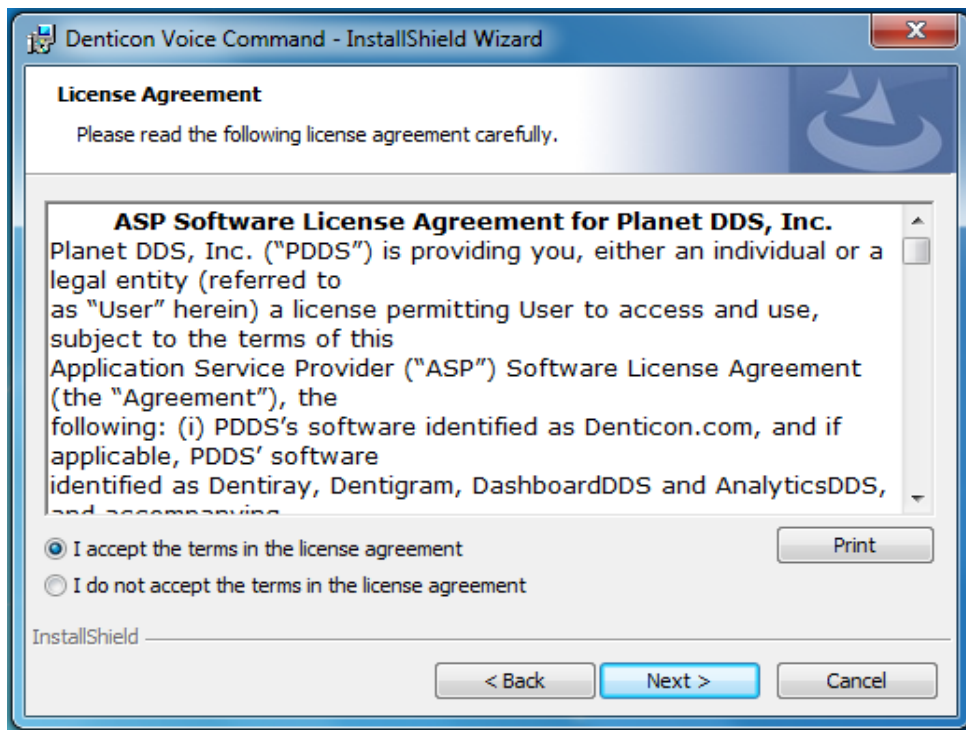

6. Click **Install**.

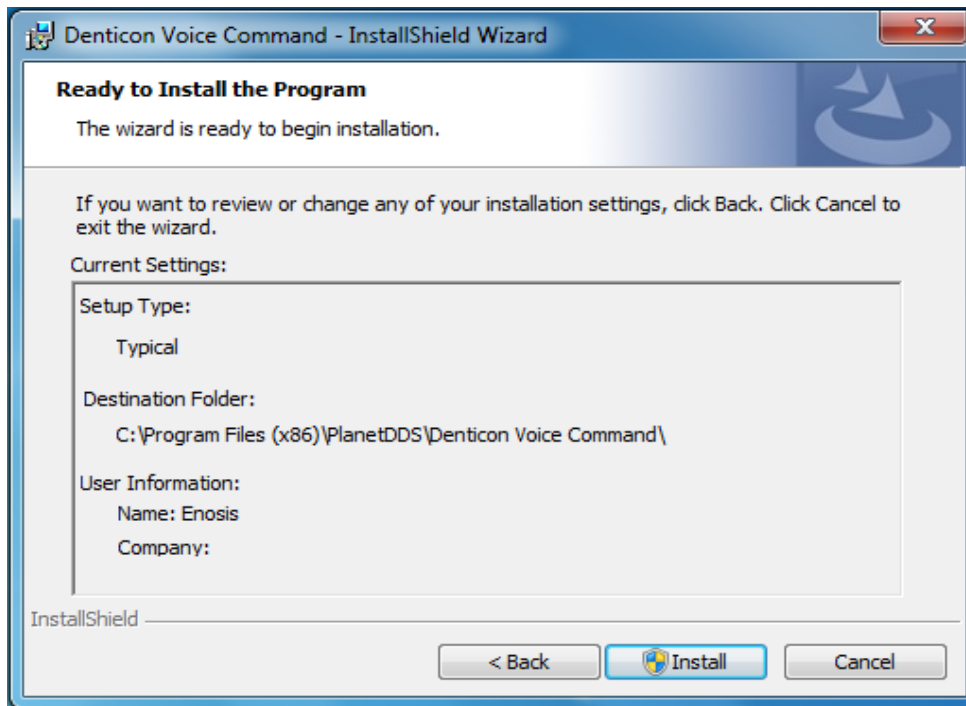

7. Click **Finish** to complete the installation

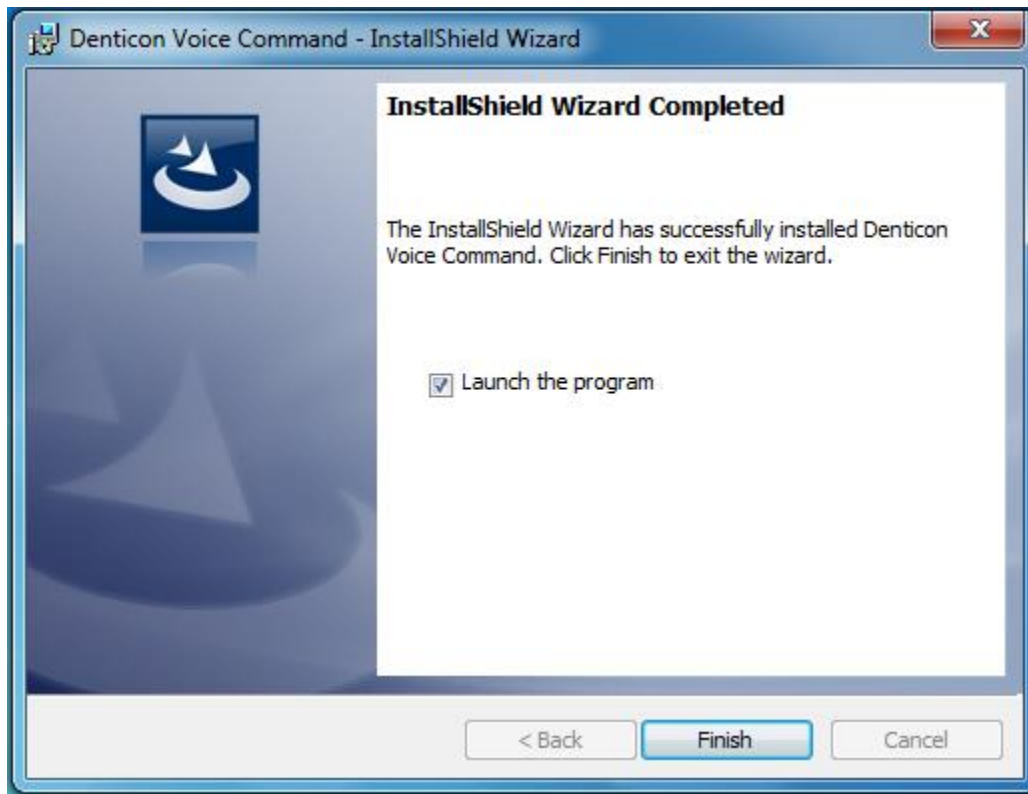


8. Now the **Denticon Voice Command** app is installed in the machine.


## Requirements

### Internet Explorer

1. Click on **Settings** -> **Internet Options**.

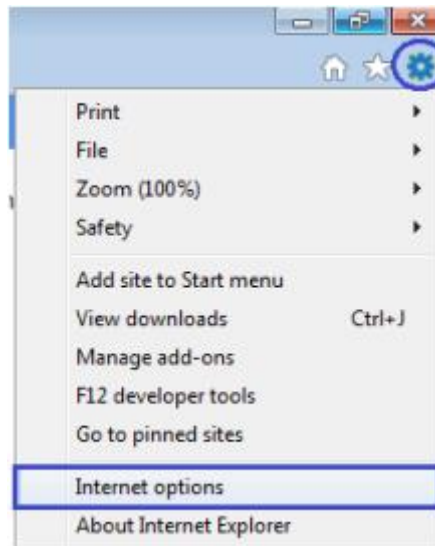

2. Then Click on **Security** tab -> **Local Intranet** tab -> **Sites** button

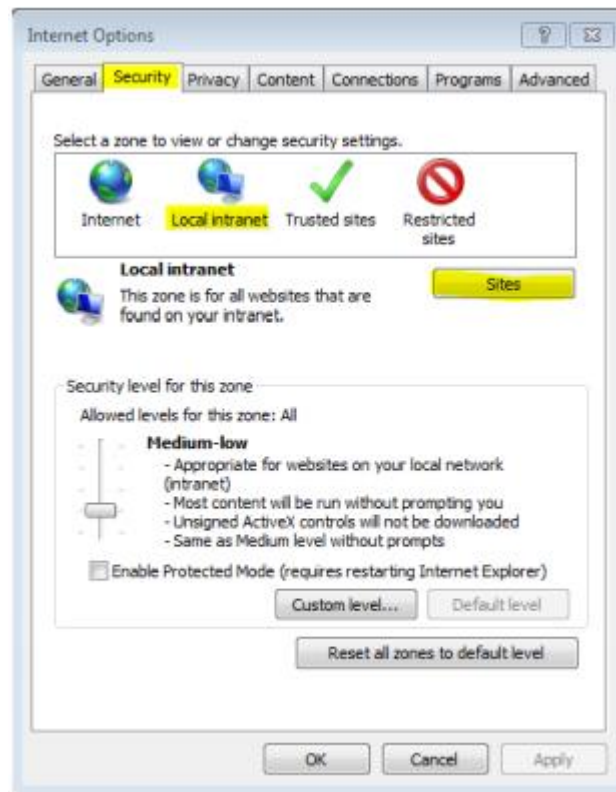


3. Now, uncheck all the boxes and press **OK**.




## Firefox

1. Go to **Menu** and click **Options**.

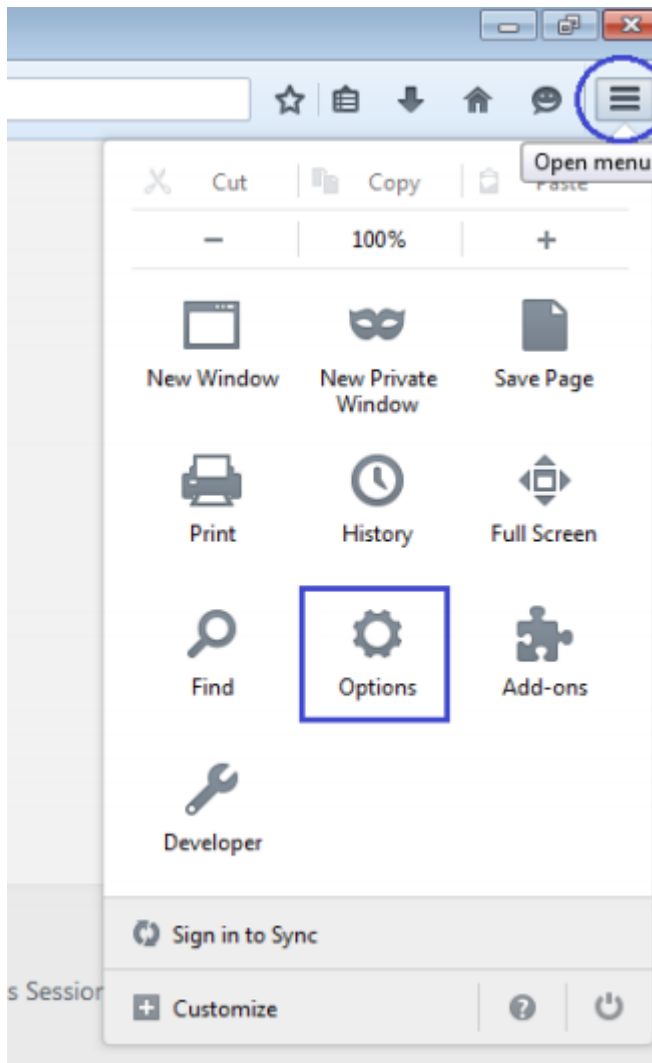





- From the Options menu, go to **Advanced** tab. Then click of **Certificates** tab and click on **View Certificates** button.

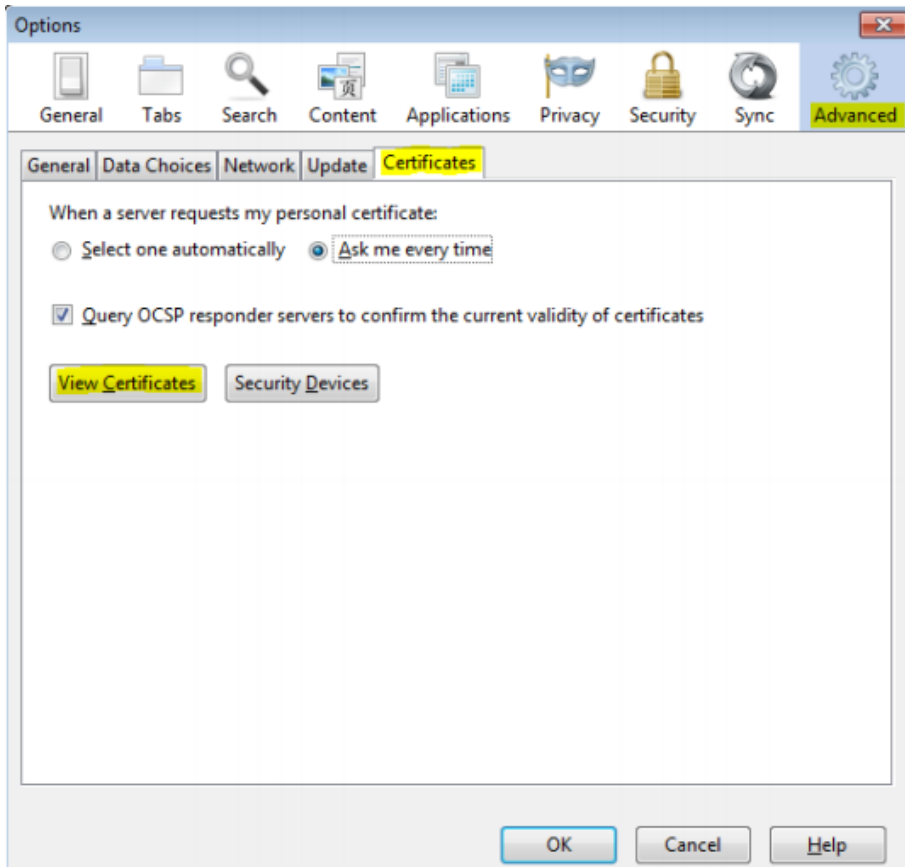


- Now, click on **Servers** tab and click on **Add Exception** button.

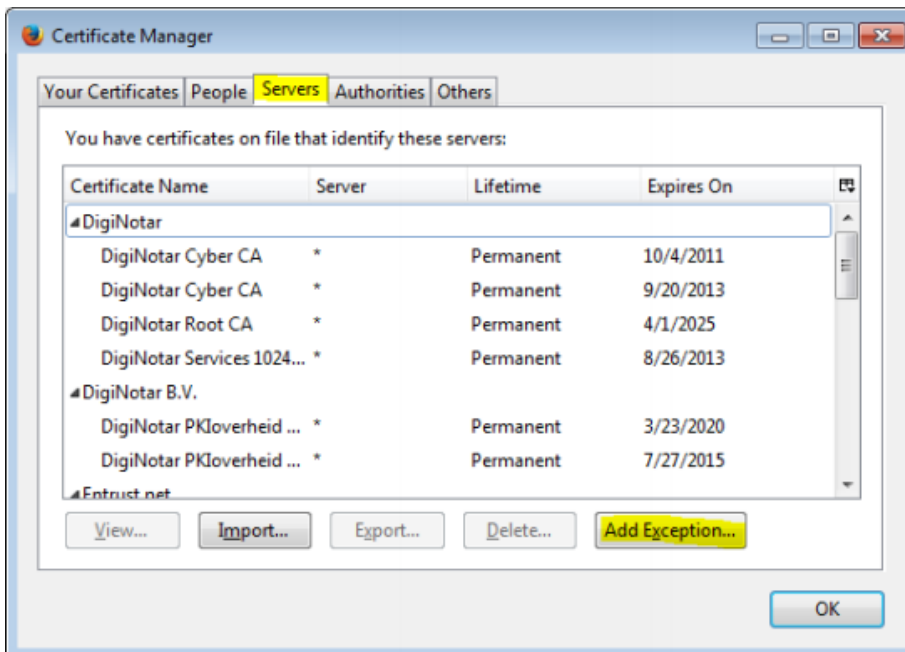

- Now, write **https://localhost:1894** in Location textbox and click on **Get Certificate**.

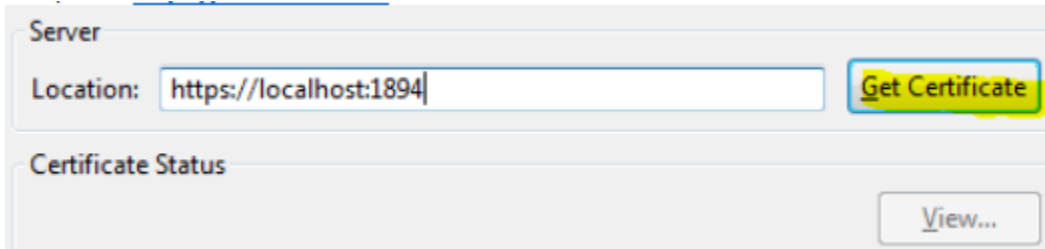

- Then, check **Permanently store this exception** and click **Confirm Security Exception**.

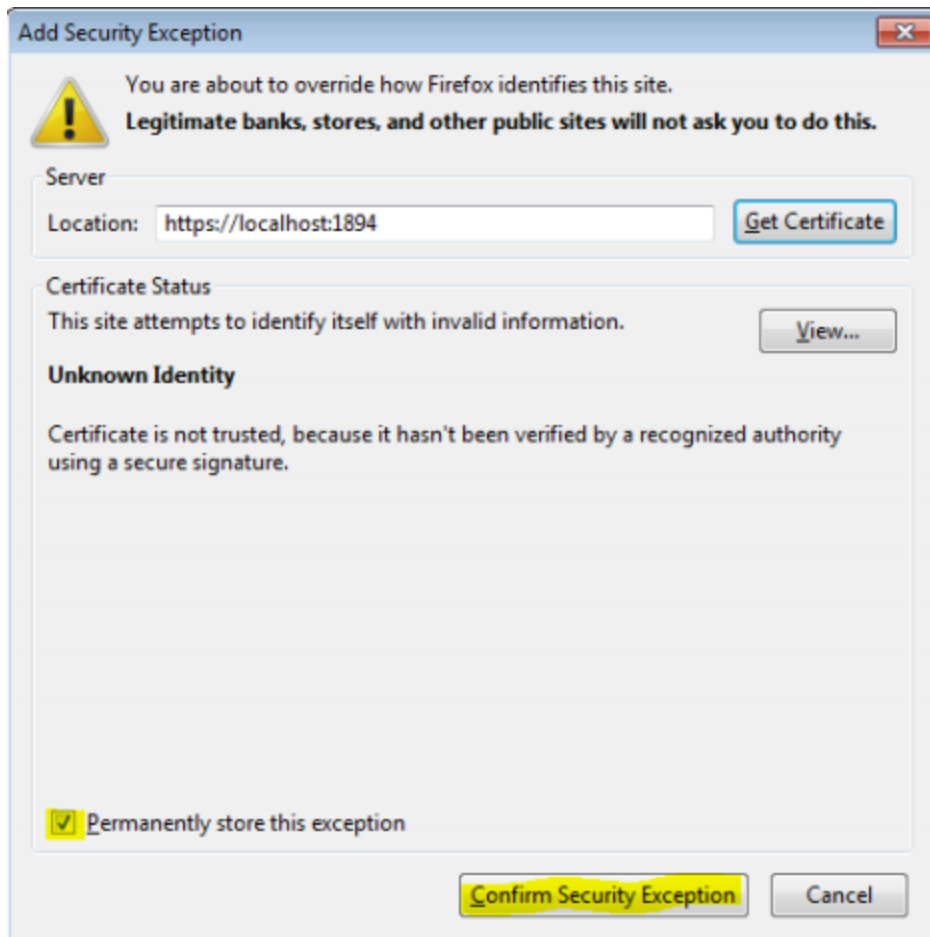


## Chrome

- No setting change required for Chrome.


# User Guide

## How to Configure

1. Login to [Denticon](#).
2. Navigate to **Setup -> Charting -> Perio Setup Templates**.
3. Select the **Perio Template Name**.
4. Click **Edit** button.

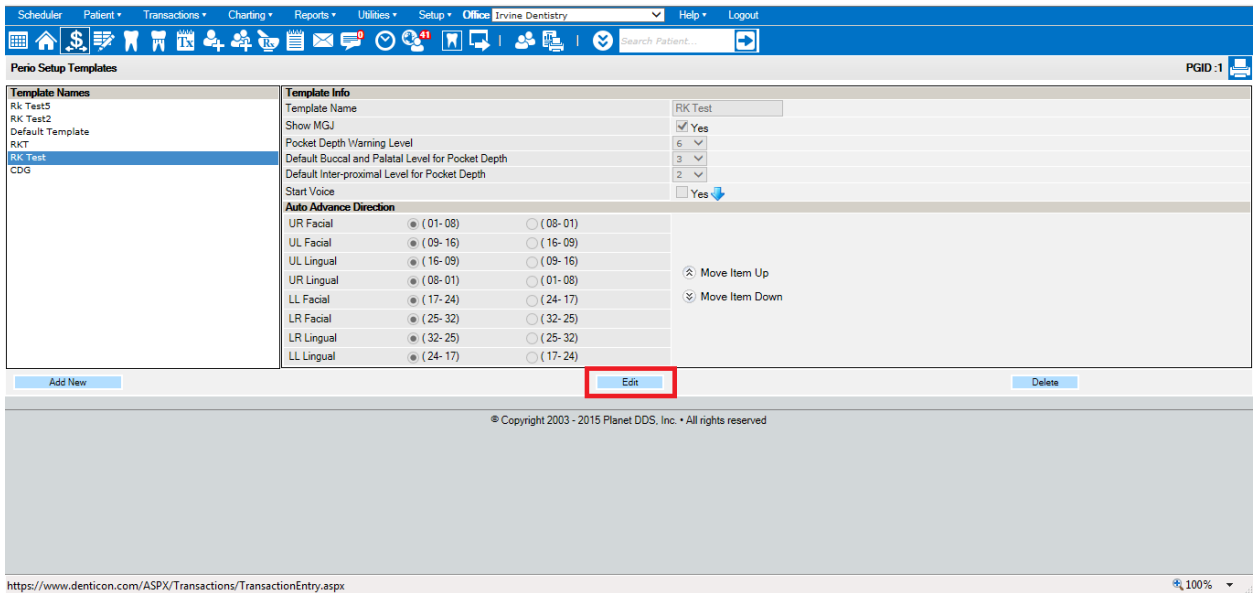


5. Check **Start Voice** and Click **Save/Update**

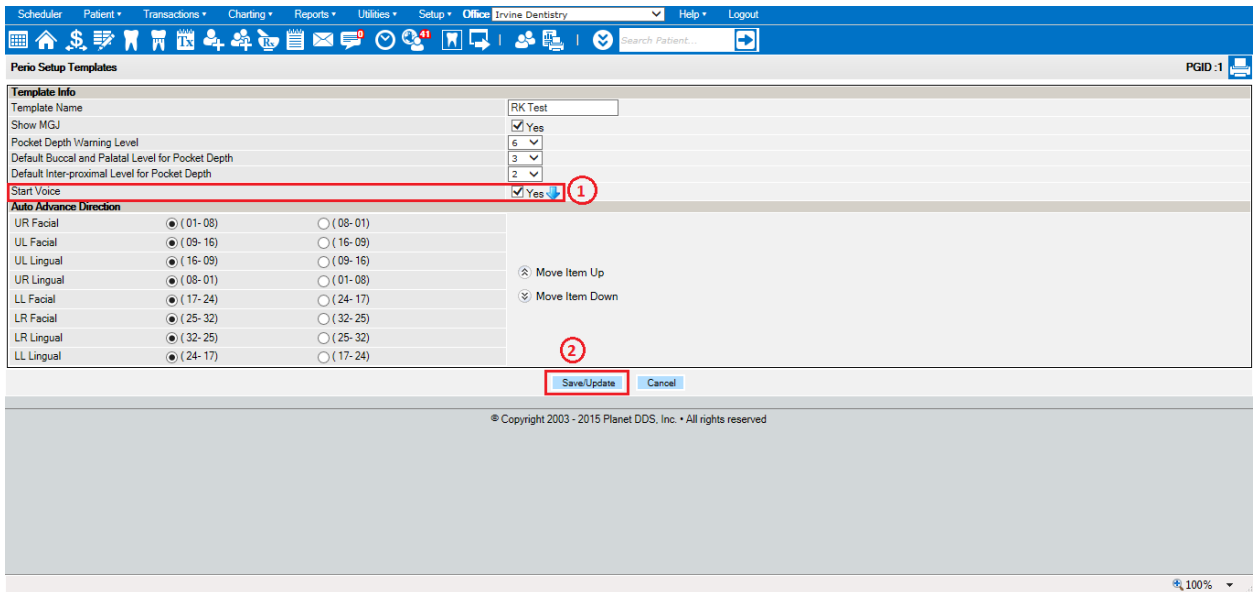

6. If **Start Voice** Checkbox is not available please contact **Denticon Support Team**.

## How to Run

1. Start **Denticon Voice Command** app from notification area

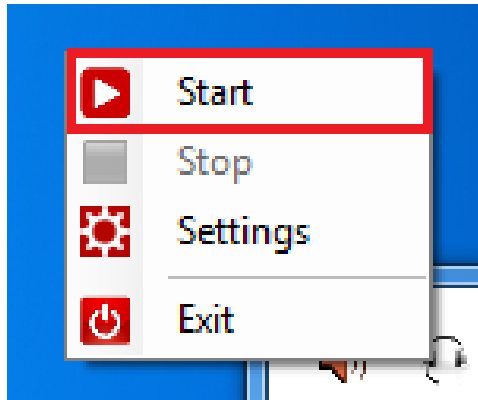

2. Navigate to **Charting-> Perio Chart (new)**.
3. Click **Start Voice**.

The screenshot shows the Denticon software interface. The top menu bar includes Scheduler, Patient, Transactions, Charting, Reports, Utilities, Setup, Office, Drs. Lee, Gutierrez, Khot And K., Help, and Logout. The patient chart area displays patient information for Alice Jacobs2, including ID, age, sex, and dates. The main area shows a Perio Chart with a grid of tooth diagrams. The right-hand panel contains a 'Start Voice' button, which is highlighted in red. Other buttons in the panel include 'New Exam', 'Exam Details', 'Compare By Dates', 'Pocket Depth', 'Set Defaults', and 'Disconnected...'. The 'Set Defaults' section includes a table with columns for 'Prev' and 'Next' values.

	Prev	Next
FGM	0	1
MGJ	2	3
Bleeding	4	5
	6	7
Suppuration	8	9
	10	11
Furcation	12	13
	14	15

© Copyright 2003 - 2015 Planet DDS, Inc. • All rights reserved



## How to Use

After the connection is made, Perio Voice should activate. Following voice commands (highlighted in green) are available in Denticon: (e.g. you can say “Pocket Depth 2 3 2”)

### ***New Exam***

For creating a new exam.

### ***Pocket Depth***

For selecting Pocket Depth. For probing, command **0 to 15** for corresponding Pocket Depth level.

### ***FGM***

For selecting FGM. For probing, command **+1 to +5, 0, 1 to 8** for corresponding FGM level.

### ***MGJ***

For selecting MGJ. For probing, command **0 to 15** for corresponding MGJ level.

### ***Bleeding***

For selecting Bleeding. Command **On** for apply bleeding and **Off** for remove Bleeding.

### ***Suppuration***

For selecting Suppuration. Command **On** for apply bleeding and **Off** for remove Suppuration.

### ***Furcation***

For selecting Furcation. Command **1 to 4** for probing corresponding Furcation Class. Command **0** for reset.

### ***Mobility***

For selecting Furcation. Command **0.5(zero point five/ point five), 1, 1.5(one point five), 2, 2.5(two point five)** for probing corresponding Mobility. Command **0** for reset.

### ***Next***

For selecting next tooth block according to the direction specified in the template.

### ***Previous***

For selecting previous tooth block according to the direction specified in the template.

### ***Jump to #Tooth No***

For jumping to the first block of the upper arch of a Specific Tooth.

### ***Go lower***

For selecting the corresponding tooth block of the selected Tooth in lower arch.


### ***Go upper***

For selecting the corresponding tooth block of the selected Tooth in upper arch.


## Voice Feedback

There's an option to enable voice feedback feature. Enabling Voice Feedback option user can hear numbers repeated back. For example for Voice Command Select **Pocket Depth**> Speak "1", Then user will hear a voice saying "1".

Enable this feature by going to Settings Menu.

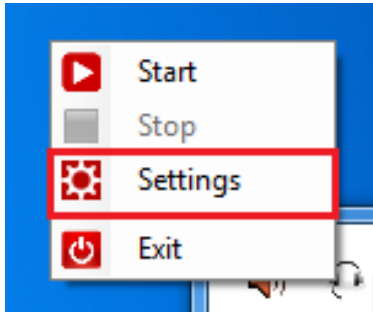

A pop up will appear, now tick the Enable Voice Feedback checkbox. And press the save button.

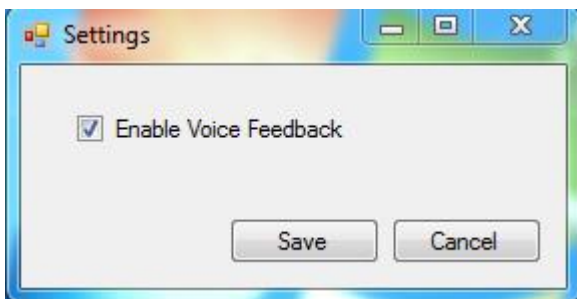


To Disable Voice Feedback option, simple uncheck the checkbox and save.

**NB: Voice Feedback only repeats back the number commands like 1, 2, 3, 4, +1, +2, +3 etc.**


## New updated version install Instructions

Follow the instructions to install new updated version.

1. If there is a new update version, download the updated new installer from Help & Support.
2. If already installed a previous version of voice command, Stop and Exit the application.

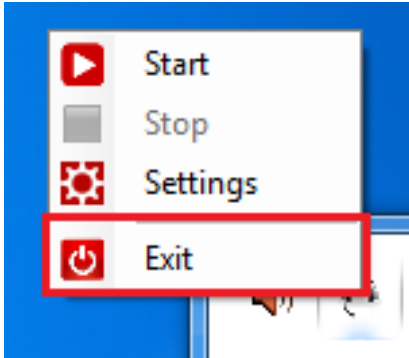

3. Uninstall the previous version from Control Panel.

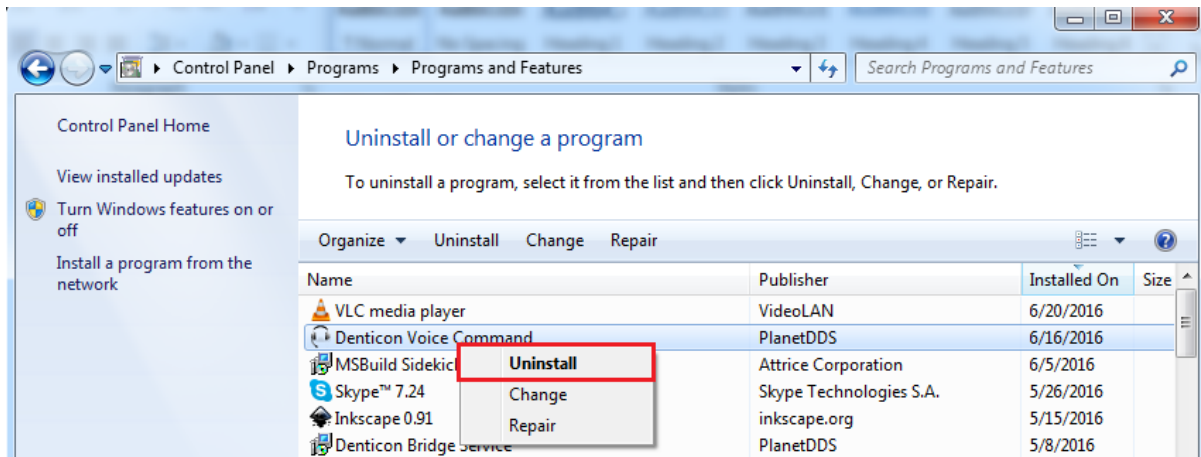

4. Restart the computer.
5. Install the latest downloaded installer.